# beurer FC 45

**ENGLISH**

**Facial Brush**

Instructions for use.....2

**ESPAÑOL**

**Cepillo facial**

Manual de instrucciones.....12

**FRANÇAIS**

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Mode d'emploi .....21

Distributed by/Distribuido por/

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Beurer North America LP

900 N Federal Highway, Suite 300

Hallandale Beach, FL 33009

[www.beurer.com](http://www.beurer.com)

**Questions or comments?**

Call toll free 1-800-536-0366 or contact [info@beurer.com](mailto:info@beurer.com)

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Llame gratis al 1-800-536-0366 o escriba a [info@beurer.com](mailto:info@beurer.com)

**Questions ou commentaires?**

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**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE  
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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# ENGLISH

## CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS.....	3	6. To Use .....	7
2. Introduction .....	5	7. Care and Maintenance .....	7
3. Package Contents .....	6	8. Specifications .....	8
4. Parts and Controls.....	6	9. FCC Part 15, Subpart B.....	8
5. Before First Use.....	6	10. Warranty and Return Policy .....	9

## SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

### 1. IMPORTANT SAFETY INSTRUCTIONS

#### Signs and symbols

The following signs appear in the Safety Section (page 3) and in this manual on page 4.

**ENSURE YOU READ AND UNDERSTAND THIS ENTIRE MANUAL BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID RISK OF INJURY.**

**WARNING**

indicates a hazardous situation which, could result in death or serious injury.

Suitable for use in bath or shower

## **WARNING:**

- The device is for adult use only and only intended for domestic/private use, not for commercial use.
- The device is only to be used on the face (avoiding eye and mouth area). Do not use the device on animals.
- The device must be used only in accordance with these instructions. The manufacturer is not liable for damage resulting from improper or careless use.
- This product is not a toy and it is not intended for use by children.
- Do not drop or hit.
- For use on face only (except for eye and mouth area).
- Do NOT use the device:
  - On cracked skin
  - On open wounds
  - In the eye area or on the mouth
  - If you have a skin disease (e.g. herpes, vascular disease or skin growths)
  - If you have skin cancer
  - If you have a skin irritation
  - Immediately after sunbathing or if you have sunburn
  - If you are taking steroidal medication
  - If you are taking blood thinners
  - If you have multiple sclerosis, diabetes or an autoimmune disorder

## Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for an extended period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge non rechargeable batteries; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose of or recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.


## 2. Introduction

The FC 45 facial brush provides gentle, deep-pore facial cleansing.

### 3. Package Contents

- 1 x Facial brush
- 1 x Replacement brush head
- 2 x AAA 1.5 V Batteries
- 1 x Instruction manual

### 4. Parts and Controls
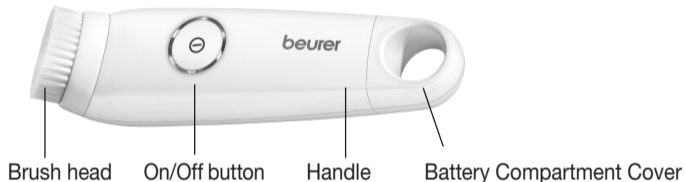
### 5. Before First Use

1. Pull the Battery Compartment Cover away from the Handle to remove. Insert the two AAA batteries (included) into the Handle according to the polarity markings inside the battery compartment.
2. Reattach the Battery Compartment Cover, ensuring it is inserted fully. If it does not fully insert, turn it over and try again.

## **6. To Use**

1. Moisten the Brush Head with water. Never use the Brush Head while it is dry as this can irritate the skin.
2. Apply cleanser or product of your choice on your face (not included).
3. The device has both a fast and a slow speed setting. Press down on the On/Off Button for 2 seconds to activate low speed. Press down on the On/Off Button again to switch to high speed.
4. Place the brush against your face and use a gentle, circular motion to remove residue, dirt and dead skin cells. To avoid irritation, do not press the brush too firmly to the skin. Treat each section of your face for a maximum of 20 seconds.
5. To shut the device off, press and hold down on the On/Off button for 2 seconds until the motor stops.
6. Cleanse your face thoroughly with water to remove any residue from your skin.
7. Pat the skin dry with a soft, dry towel and apply a moisturizing cream of your choice.
8. Recommended usage time is one minute per session. The device will shut off automatically after one minute.

## **7. Care and Maintenance**

- Before cleaning, pull the brush head off the handle. Clean the brush thoroughly with water after each use.
- Protect the product from knocks, dust, chemicals, extreme temperatures and nearby sources of heat (ovens, heaters).
- Clean the Handle with a soft, slightly moistened cloth. You can also moisten the cloth with mild soapy water. Do not use any abrasive cleaning products.

## 8. Specifications

Dimensions	6.3 x 1.6 x 1.4 in (16 x 4 x 3.5 cm)
Weight	approx. 3.3 oz (94 g) (without batteries)

## 9. FCC Part 15, Subpart B

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## 10. Warranty and Return Policy

### Limited Lifetime Warranty For Original Purchaser

This beurer Facial Brush, Model FC45, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the beurer Facial Brush, Model FC45, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitution in lieu of repair.

For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of

Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is operated in the country in which it was purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

**THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein. This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any

products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

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