beurer

IPL 7500

ENGLISH	IPL 7500 Permanent hair reduction Instruction manual	3
ESPAÑOL	IPL 7500 Reducción permanente del vello Manual de instrucciones	21
FRANÇAIS	IPL 7500 Réduction permanente de la pilosité Mode d'emploi	41

Distributed by/Distribuido por/

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Beurer North America LP
900 N Federal Highway, Suite 300
Hallandale Beach, FL 33009
www.beurer.com

Questions or comments?

Call toll free 1-800-536-0366 or contact info@beurer.com

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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ENGLISH

Contents

1. IMPORTANT SAFETY NOTES	3	10. FCC Statement	17
2. Important information	5	11. Warranty	18
3. Introduction	6		
How IPL (intense pulsed light) technology works	6		
Long-term hair reduction.....	6		
4. Parts and Controls	7		
Who can use the device.....	9		
Who cannot use the device.....	9		
Possible side effects	10		
Benefits	11		
5. To Set Up	12		
Preparing Skin for Treatment.....	12		
Preliminary Test	13		
Treating the Skin.....	13		
Recommended Treatment Schedule.....	14		
What to expect when using the device	15		
After Treatment.....	16		
6. Care and Maintenance	16		
7. Troubleshooting	16		
8. Specifications	17		
9. Storage	17		

1. IMPORTANT SAFETY NOTES

Signs and symbols

The following signs appear in the Safety Section (page 4) and in this manual on pages 7, 10, 13 and 14.

WARNING

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE

NOTICE addresses practices not related to personal injury, such as product and/or property damage.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

WARNING

To reduce the risk of fire, electric shock, or serious personal injury:

- Do not open or repair the device as it may expose you to dangerous electrical components which may cause serious bodily damage.
- Do not place, store or use the device near bathtubs, showers, basins or other vessels containing water.
- Do not leave the device unattended while plugged in.
- Never look directly at the light coming from the device applicator and lamp.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin one hour before using the device.
- Do not use the device if you have an active implant such as a pacemaker, insulin pumps, or similar devices.

- Do not use the device on tanned skin or after sun exposure. Avoid exposure to the sun for four weeks before using this device.
- Do not use the device on naturally dark skin (see guidance at the back cover and high Fitzpatrick level in the “Who should use the device” section).
- Do not use the device if you have dark brown or black spots such as freckles, birth marks, around the nipples, moles or blisters on the area to be treated.
- Do not use the device on mucous membranes such as nostrils, lips, eyelids, ears, genital areas or anus.
- Remove all piercings or metal items such as earrings or jewelry before treatment.
- Do not use the device on tattoos or permanent make-up.

CAUTION

To reduce the risk of personal injury or product/property damage:

- Use the device only according to the instructions for use.
- Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk group for

this response is females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.

- This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities.
- Do not use the device if:
- The device is not working properly, appears damaged or broken (UV and IR emitted from this device)
- The cord or plug are damaged
- You see or smell smoke
- The outer shell is damaged or loose
- Lamp or Treatment Window is damaged/missing
- Optical skin tone is cracked, or broken
- Becomes damp or wet.
- Avoid more than 15 minutes of sun exposure for at least 2 weeks post treatment. If during these 2 weeks your treated skin area is exposed to sunlight, make sure to apply SPF 30 or higher sunscreens.
- Do not treat the same area of skin more than once per treatment as it may cause skin irritation.
- **Use the device only with its original power supply. For your safety, never use DC adapters which are not approved by the manufacturer.**

- **This device is intended for one individual's use only. To help prevent transmitting diseases between users, do not let others use the device.**
- **If you have undergone any other aesthetic procedure or any other uncertainty regarding potential side effects, consult your physician before using the IPL 7500.**

2. Important information

⚠ Signs and symbols

The following symbols are used in these instructions for use, on the packaging and on the type plate for the device and accessories:

	Caution
	Note Note on important information
	Follow instructions for use

	Type BF applied part
	Device in protection class 2
IP22	Device protected against foreign objects \geq 12.5 mm and against falling drops of water
	Permissible storage temperature and humidity
	Permissible operating temperature and humidity
	Keep away from eyes
	Do not use the IPL 7500 device on tanned skin or after sun exposure
	Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).
	Keep dry / Protect from moisture / Keep away from showers, basins or other vessels containing water.

	Use device only with its original power supply. For your safety, never use DC adapters which are not approved by the manufacturer
--	---

3. Introduction

The device is an over-the-counter device designed to remove unwanted hair and permanently reduce regrowth as measured at 6, 9, and 12 months after treatment ends.

How IPL (intense pulsed light) technology works

The device is based on the principle of selective photothermolysis. IPL uses a combination of Intense Pulsed Light (IPL) technology and unique skin condition sensors. The light pulse is absorbed by the melanin in the hair shaft, causing local thermal damage which results in a reduction in hair growth.

Long-term hair reduction

The hair growth cycle consists of three phases: anagen, catagen, and telogen. Anagen is the growth phase, while

catagen and telogen are resting phases. Permanent hair reduction can only occur during anagen.


Skin types are depicted on the Fitzpatrick scale below. Darker skin contains more melanin, which competes with melanin in the targeted hair for energy absorption. When more light energy is absorbed into the skin, there is a greater risk of adverse effects and epidermis damage. This device is only recommended for use on skin types I through IV.

WARNING

- Read this user guide completely before starting treatment with the device.

4. Parts and Controls

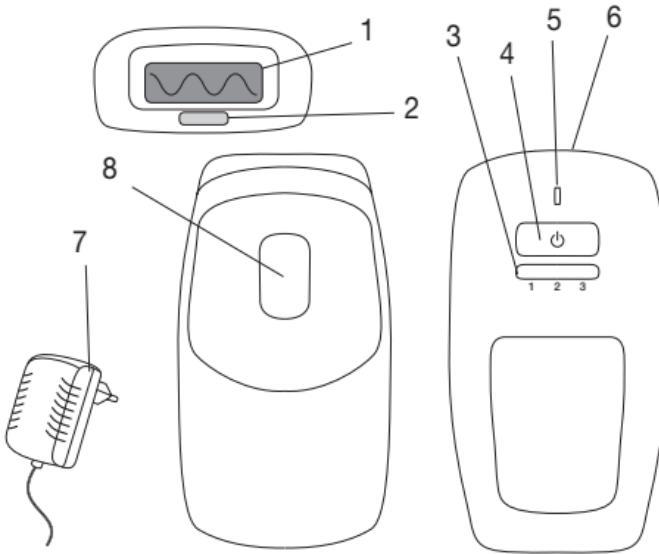1	Light area/treatment window/lamp
2	Skin color sensor/skin contact sensor
3	Intensity display (3 levels)
4	ON/ OFF power/ energy selection button
5	Indicator lamp (indicating lack of contact/not allowed skin tone/device not usable)

6	Socket for device plug
7	Mains adapter
8	Trigger button

ON/ OFF power/ energy selection button

- Press the ON/ OFF power/ energy selection button (after the device is connected to the power supply) to activate device at Level 1. Each press will increase the energy level from 1 through 3.
- Press and hold the ON/ OFF power/ energy selection button for three seconds to shut the device off. The device will also shut off by itself if not touched for 180 seconds.

The Trigger button (flash pulse activation) activates the flash lamp as long as the treatment window is in full contact with the skin and the skin tone is in the allowed range. Press the trigger button once to emit a single pulse; press and hold the trigger button for continuous pulses.

Indicator	Mode	Description
No light	OFF	Device not working (not connected to the power supply)
Blinking green light	Standby	Device is in standby mode (when first powering on or when no buttons are pressed for 3 minutes)
Green energy level LEDs (1 -3)	Energy level 1-3	The three energy levels are indicated by three green LEDs.
Steady green light	ON	Pulse is ready
Short orange light flash	Activation (pulse) blocked	The device does not deliver a pulse due to the lack of full skin contact.
Short red light flash	Activation (pulse) blocked	The device does not deliver a pulse due to non-allowed skin tone.

Blinking red light	Close to end of life	Pulse light is nearing or reached the end of its service life.
"Steady" red light	General Error	Device has malfunctioned.

Who can use the device

The device is intended for women and men who desire to remove unwanted hair.

This device can be used on all areas of normal skin on the body including legs, underarms, bikini line, back, and arms. This device should not be used on sensitive areas such as genitals, lips, mouth or near the eyes.

This device can be used on skin tone types 1 to 4 as depicted in the Fitzpatrick skin chart below. Do not use on dark skin such as types 5 and 6.

The device works best on dark hair, which contains more melanin. Lighter hair will respond, but typically will require more treatment sessions. White, gray or blond hair shows a more limited response, but some users have noted hair reduction after multiple treatment sessions.

Never use the device around or near your eyes.

Who cannot use the device

To prevent the possibility of serious personal injury, do not use this device in the following instances, since the safety of this device in these cases has not been fully established.

- If you are under the age of 18.
- If you are pregnant, nursing or might be pregnant.
- If you have skin cancer, on areas at risk of malignancy, or if you have received radiation therapy or chemotherapy within the past three months.
- If you have epilepsy.
- On damaged or broken skin or if you have a history of herpes or psoriasis outbreaks in the treatment area.
- If you suffer from keloid scar formation (growth of extra scar tissue where the skin has healed after an injury).

- If you have been taking tetracycline antibiotics and medication with the following active ingredients within the last two weeks: Retinoid acid, Retin-A, or drugs containing Vitamin A such as Accutane.
- If you are photosensitive or have other diseases related to light sensitivity.
- If you suffer from severe or decompensated metabolic conditions such as diabetes since it may delay healing and infection in the case of an adverse effect.
- If you have had a skin peel or other skin resurfacing procedure within the last eight weeks.

CAUTION

- If you have undergone any other aesthetic procedures or any uncertainty regarding potential side effects, consult your physician before use.

Possible side effects

If used according to these instructions, side effects from this device are rare. Any aesthetic procedure can cause some degree of side effects, however. It is important to understand the risks and complications that could occur with a pulsed light hair removal system.

Degree of Side Effects	Side Effects
Minor	Skin Irritation Skin Redness Skin Sensitivity
Rare	Skin Burns Scarring Pigment Changes Excessive Redness and Swelling
Extremely Rare	Infection Bruising

- **Skin Irritation**

The use of Pulsed Light Hair Reduction System is generally very well-tolerated, but many users do feel a sensation of warmth and a mild stinging sensation during treatment or a few minutes after. Anything beyond minor skin irritation means that either the energy level is too high or your skin is unable to tolerate the hair reduction process.

- **Skin Redness**

Redness on the treated area may appear after treatment or within 24 hours. Consult your physician in case the redness does not disappear after 2 to 3 days.

- **Skin Sensitivity**

You may feel that skin in the treated area is more sensitive or you may encounter dryness or flaking of the skin.

- **Skin Burns**

Rarely, burns on the treated area can occur following the use of the hair removal reduction system. The burns can require a few weeks to heal and, extremely rarely, may leave a visible scar.

- **Scarring**

In rare cases, permanent scarring may occur. Usually when scarring occurs it is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

- **Pigment Changes**

Use of the device may cause temporary skin hyperpigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) to the surrounding skin.

The risk is higher for people with darker skin tones. Usually discolouration or changes to skin pigment are temporary and permanent hyperpigmentation or hypopigmentation rarely occur.

- **Excessive Redness and Swelling**

Very rarely, the treated area skin may become red and swollen. This side effect is most common in sensitive areas of the body. Redness and swelling should decrease after 2 to 7 days and you should treat the area by applying ice. Avoid exposure to sun.

- **Infection**

In rare cases, infection of the skin could occur after a burn caused by the treatment.

- **Bruising**

Blue-purple bruising is exceedingly rare and could last 5 to 10 days. As bruising fades, there could be a permanent rusty brown discolouration of the skin (hyperpigmentation).

Benefits

This device features a built-in skin color sensor that continually evaluates your skin tone and allows treatment only on suitable areas. This safety feature also disables treatment if the skin in the treatment area is too dark or tanned,

even if your general skin type is lighter. As an extra precaution, this device only emits a light pulse when it is in contact with your skin.

Response to the device hair reduction differs among individuals, depending on biological factors and how closely treatment instructions are followed. Based on the treatment regime in these instructions, some results can be seen in 3-4 weeks.

The hair will continue to grow up to two weeks after a treatment session. Since hair grows in three different stages, the hairs will naturally and gradually fall out depending on the hair growth stage in which the treatment was done (only hairs in the active growth stage will be affected). Some hair may not be affected by the device because of the different hair growth stages. Continue treatment over a period of time to ensure hair is treated in the effective stage.

For permanent reduction in hair regrowth, defined as a long-term, stable reduction in the number of hair regrowing when measured at 6, 9, and 12 months you should follow treatment regimen.

5. To Set Up

1. Remove all components from the box (handheld device and power supply).

2. Before beginning the treatment, make sure the device is clean, dry, and intact.
3. Plug the power supply end into an electrical outlet and the cord end into the device socket. The device is now ready to be used.

Preparing Skin for Treatment

1. Before treatment, make sure your skin is shaved, clean, and dry (free of any residue such as powders, creams, makeup, antiperspirants, or deodorants).
2. Between treatments, you can shave but do not wax or pluck treated hair.
3. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin for at least one hour before treatment.
4. Ensure the device is clean and then perform a preliminary test on a small patch of skin in the intended treatment area 48 hours prior to first full use. This test will help determine the proper energy level setting for your skin type (see the “Preliminary Test” section below). The device is now ready to be used.
- For a prolonged use, device cleaning is a must (see chapter 6 “Care and maintenance”).

- Remember: The device comes with a built-in skin-tone sensor, enabling the application only on a suitable skin range. This unique safety feature will not let you treat areas where your skin is too dark or too tanned, even if your general skin type is lighter.
- The device is suitable for men. However male body hair, especially those on the chest will require more hair reduction sessions.

Preliminary Test

A preliminary test will determine your skin's reaction to the device treatment and help you set the optimal energy level for each treatment area. Perform a preliminary test on each area you intend to treat 48 hours prior to the full treatment.

CAUTION

- **Do not treat full areas with energy levels that you have not tested before! Do not treat the tested area for at least 48 hours after the preliminary test is completed.**
1. Select an area of skin to test and place the treatment window in full contact with the skin. Press the On/Off/power selection button once to select Energy Level 1 and then press the Trigger button once to emit a single pulse.

2. If you feel comfortable with the energy level, press the On/Off/power selection button again to increase to Energy Level 2, move the treatment window of the device to another test area of skin and press the trigger button again.
3. If you still feel comfortable, increase the Energy Level to 3.
4. Stop the test if you do not feel comfortable, as increasing the energy level increases the risk of adverse effects.
5. If you see adverse effects such as severe redness or blistering, stop the test immediately.
6. Wait 48 hours and examine the tested area. If the skin in the test area seems normal (no response or mild redness observed), you can proceed with a full treatment in the highest comfortable energy level. If you see excess redness, swelling, or blistering even in the low-energy level settings, do not treat that area.

Treating the Skin

IMPORTANT:

Conduct a preliminary test before performing a full treatment. See the "Preliminary Test" section above.

1. Plug the power cord into the device; the green power LED will flash indicating the device is in standby mode. Press the On/Off/power selection button to turn on the

device; one green Level LED will illuminate indicating the device is set to Energy Level 1. To increase the Energy Level, press the On/Off/power selection button repeatedly to select Energy Level 2 or 3.

Select the energy level according to your convenience, as determined in the preliminary test for that treatment area.

2. Position the treatment window onto the desired area of skin, making sure it makes full contact. If the treatment window is not properly positioned, the Power LED will briefly change to orange. If your skin tone is too dark/tanned, the Power LED will briefly change to red. In either case, the device will not operate.
3. Press the trigger button for a single pulse or press and hold for continuous operation. After the first pulse, the internal fan in the device may switch on.
4. When pressing the trigger button you will see a bright flash, hear a popping noise, and may feel a mild stinging and/or heating sensation.
5. Keep the trigger button pressed and move the device to a different area after each pulse while keeping it in contact with the skin. Continue until the entire treatment area has been covered.
6. After treatment, shut the device off by pressing and holding the On/Off/power selection button for three seconds;

the orange LED will illuminate while the device is in shutdown mode, followed by blinking green which indicates the device is in standby mode. Unplug the device after each treatment session.

CAUTION

- **Do not overlap pulses in the same spot. Repeated pulses on the same spot may cause burns.**
- **Do not treat the same area of skin more than once per treatment as it may cause skin irritation.**
- **The device will not trigger a pulse if:**
 - a. There is no full contact of the treatment window with the body skin.
 - b. Your skin tone is too dark or tanned.

Recommended Treatment Schedule

Body Hair Color	Light to Medium Brown	Dark Brown to Black
Fitzpatrick Skin Tone	1-4	1-4

Recommended Energy Level	1-3	1-3
Treatment Area	Legs	Legs
	Arms	Arms
	Underarms	Underarms
	Bikini Lines	Bikini Lines
Average Number of Treatments	10-12	8-10
Treatment Interval	The first four treatments should be once every two weeks, then four treatments, once every four weeks.	

NOTES:

- Hair reduction effectiveness differs among individuals, depending on biological factors.
- It may take up to 3-4 weeks before full results can be seen. Hair will continue to grow up to two weeks after a hair reduction session. Since hair grows in three different stages, the hairs will naturally and gradually fall out depending on the hair growth stage in which

the treatment was performed (only hair in the active growth phase will be affected). Some hair may not be affected because of the different hair growth phases. To remove the largest amount of hair, multiple treatment sessions are required.

- The number of sessions required is different for each treatment area and for each user, as it is affected by type of skin, hair, and hormone levels.
- Maintenance may be needed from time to time if growth is still visible.

What to expect when using the device

During treatment it is normal to experience:

- A flash of light on each press of the trigger button. The light will not harm the eyes when the device is applied with full contact to the skin on the recommended treatment areas.
- A popping noise with each pulse.
- A sensation of heat radiating from the treated area after each pulse.
- A mild stinging sensation, similar to being snapped with a rubber band.
- Some mild red or pink coloration during and just after treatment.

After Treatment

- Apply a soothing lotion such as aloe vera.
- Do not use peeling and bleaching creams or products for at least 24 hours.
- Avoid direct sun exposure for at least two weeks. If recently-treated skin is exposed to sunlight in this time period, make sure to apply SPF 30 or higher sunscreen.
- Do not wax or pluck hairs from treated areas.

Remember that hair may take up to two weeks to fall out if the hair follicle is still active.

6. Care and Maintenance

- Shut off and unplug the device before cleaning.
- Clean the applicator head cover after every session by using a dry, clean cloth.
- Keep the device dry. Never place, store or use this device near bathtubs, showers, basins or other containers of water as this may cause electric shock.
- If you do not plan on using the device for a long period of time, clean it and store it away.
- If you plan to travel and take your device with you, you should store it in the original box and follow the maintenance instructions.

7. Troubleshooting

- If the device does not turn on:
 - a. Ensure the power cord is connected properly into the device and the electrical outlet.
 - b. Make sure that the electrical outlet has power.
- If the device does not trigger a pulse:
 - a. Ensure it makes full contact with the skin. The trigger button will operate only when the device is in full contact with the skin.
 - b. Check that the area being treated is large enough and sufficiently flat. The whole treatment window must touch the skin area. Pull the skin flat if needed.
 - c. If skin tone is too dark, the device will not operate and the red indicator LED will flash.
 - d. In the case of a safety issue or when the lamp has reached the end of its service life, a red LED will illuminate and then the device will shut down.
- General error:
 - a. Keep the treatment window clean. If you cannot get the area clean, contact Customer Service.
 - b. If pulses are delayed after pressing the trigger button, try to move to a cooler room. The device may pulse at a slower rate in a warm environment.

- c. In case of a general error, unplug the power cord from the electrical outlet, wait 10 seconds and then reconnect it to reset the device.
- If you still have problems, contact Customer Service.

8. Specifications

Emitted light spectrum	475nm–1.200nm max
Emitted energy flux:	5 Joules/0.16 sq in (5 Joules/cm ²)
Selectable energy levels:	3
Treatment window area:	0.47 sq in (3 cm ²)
Pulse Rate:	1 pulse every 3 seconds
Power Supply Output:	12 VAC, 1.8 Amps
Power Supply Input:	100–240 VAC, 50–60Hz
Dimensions:	2.83 x 5.2 x 1.73 in (72 x 132 x 44 mm)
Weight:	7.9 oz (224 g)

9. Storage

Storage and Transport Conditions	
Temperature	14 °F to 158 °F (-10 °C to 70 °C)
Humidity	0 to 90% RH at 131 °F (55 °C)
Pressure	non-condensing

Operational Conditions	
Temperature	50 °F to 104 °F (10 °C to 40 °C)
Humidity	30% to 70% RH
Max operating altitude	10.000 feet

10. FCC Statement

Changes or modifications to the product not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTICE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is en-

couraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer device for permanent hair reduction, Model IPL7500, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer device for permanent hair reduction, Model IPL7500, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this

written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond

the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein. This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer. This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
900 N Federal Hwy, Suite 300
Hallandale Beach, FL 33009. USA
www.beurer.com

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Table 1

For all me equipment and me systems

- The device needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided.
- Portable and mobile RF communications equipment can affect device.
- The use of accessories or replacement parts, other than those sold/specify by the manufacturer may result in increased emissions or decreased immunity of the device.
- The device should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the device should be observed to verify normal operation in the configuration in which it will be used.

Guidance and manufacture's declaration – electromagnetic emissions

IPL 7500 is intended for use in the electromagnetic environment specified below. The customer or the user of IPL 7500 should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	IPL 7500 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	IPL 7500 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	complies	

Table 2**For all me equipment and me systems****Guidance and manufacturer's declaration – electromagnetic immunity****Guidance and manufacturer's declaration – electromagnetic immunity**

IPL 7500 is intended for use in the electromagnetic environment specified below. The customer or the user of IPL 7500 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV input/output lines	± 2 kV for power supply lines ± 1 kV input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	± 1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5 % U_T (> 95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles < 5 % U_T (> 95 % dip in U_T) for 5 sec	< 5 % U_T (> 95 % dip in U_T) for 0,5 cycle 40 % U_T (60 % dip in U_T) for 5 cycles 70 % U_T (30 % dip in U_T) for 25 cycles < 5 % U_T (> 95 % dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of IPL 7500 requires continued operation during power mains interruptions, it is recommended that IPL 7500 is be powered from an uninterruptible power supply or a battery.
Power frequency (50/ 60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTICE: U_T is the a.c. mains voltage prior to application of the test level.

Table 3**For all me equipment and me systems that are not life-supporting****Guidance and manufacturer's declaration – electromagnetic immunity****Guidance and manufacturer's declaration – electromagnetic immunity**

IPL 7500 is intended for use in the electromagnetic environment specified below. The customer or the user of IPL 7500 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 V _{r m s} 150 kHz to 80 MHz	3V	<p>Portable and mobile RF communications equipment should be used no closer to any part of IPL 7500, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance:</p> $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P} \text{ 80 MHz to 800 MHz}$
Radiated RF IEC 61000-4-3	3 V / m 80 MHz to 2,5 GHz	3V/m	$d = 2,3 \sqrt{P} \text{ 800 MHz to 2,5 GHz}$ <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strength from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radios broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which IPL 7500 is used exceeds the applicable RF compliance level above, IPL 7500 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating IPL 7500.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.

Table 4

For all me equipment and me systems that are not life-supporting

Recommended separation distances between portable and mobile RF communications equipment and IPL 7500

Recommended separation distances between portable and mobile RF communications equipment and IPL 7500

IPL 7500 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of IPL 7500 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and IPL 7500 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2 \sqrt{P}$	80 MHz to 800 MHz $d = 1,2 \sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3 \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.